

The Role of Advanced Hemodynamic Monitoring in Enhanced Recovery After Surgery (ERAS) Initiatives

Using continuous hemodynamic monitoring as a tool to help maintain patient fluid balance and reduce postoperative complications.

INTRODUCTION

Patients undergoing major surgery are faced with an inherent risk of morbidity and mortality (See Table 1). These risks can increase depending on a patient's cardiovascular and hemodynamic condition and are known to contribute to a variety of postoperative complications and increased lengths of stay (LOS) in the hospital.^{1,2} Specific factors that influence a patient's LOS during postoperative rehabilitation include the need for analgesia, intravenous fluids, and lack of mobility.³ To minimize recovery time and reduce postoperative complications for a variety of high- to moderate-risk surgical patients, hospitals and surgical teams around the world have adopted a comprehensive set of perioperative practice guidelines known as Enhanced Recovery After Surgery (ERAS).

TABLE 1. POSTOPERATIVE COMPLICATION RATE²

SURGERY	MORBIDITY RATE
Esophagectomy	55%
Pelvic exenteration	45%
Pancreatectomy	35%
Colectomy	29%
Gastrectomy	29%
Liver resection	27%

ERAS guidelines consist of 22 preoperative, intraoperative, and postoperative protocols, which have been shown to lead to a reduction in complications and hospital LOS, improvements in cardiopulmonary function, earlier return of bowel function, enhanced mobilization, and earlier resumption of normal activities.^{4,5,6} These practice guidelines represent a fundamental shift in perioperative care for numerous types of surgical interventions including colorectal surgery, vascular surgery, thoracic surgery, radical cystectomy, and orthopedic cases.^{7,8,9,10,11}


One critical element of all ERAS programs is a protocol known as perioperative goal-directed therapy (PGDT), which helps ensure adequate hydration and maintain euvolemia, while avoiding hypervolemia or hypovolemia that can contribute to postoperative complications (See Figure 1).^{12,13} To achieve optimal fluid balance for the surgical patient, PGDT relies on continuous monitoring of a variety of hemodynamic targets, which can all be derived from minimally invasive cardiac output (CO) monitoring technologies.

PERIOPERATIVE GOAL-DIRECTED THERAPY

PGDT is an integral element of ERAS and has been shown in multiple single-center studies, quality improvement studies, and published meta analyses to significantly improve patient outcomes.^{14,15,16,17} By improving cardiovascular function and balancing fluid intake, PGDT helps clinicians maintain adequate oxygen supply perioperatively. As part of ERAS guidelines, PGDT helps decrease nausea, vomiting, and the incidence of ileus, or intestinal obstruction, while allowing patients to take solid food earlier, become more alert, and start walking sooner after surgery, ultimately reducing hospital LOS.^{12,13}

One critical element of all ERAS programs is a protocol known as perioperative goal-directed therapy (PGDT), which helps avoid hypervolemia or hypovolemia that can contribute to postoperative complications.^{12,13}

FIGURE 1. COMPLICATIONS ASSOCIATED WITH SUBOPTIMAL PERIOPERATIVE FLUID MANAGEMENT¹³

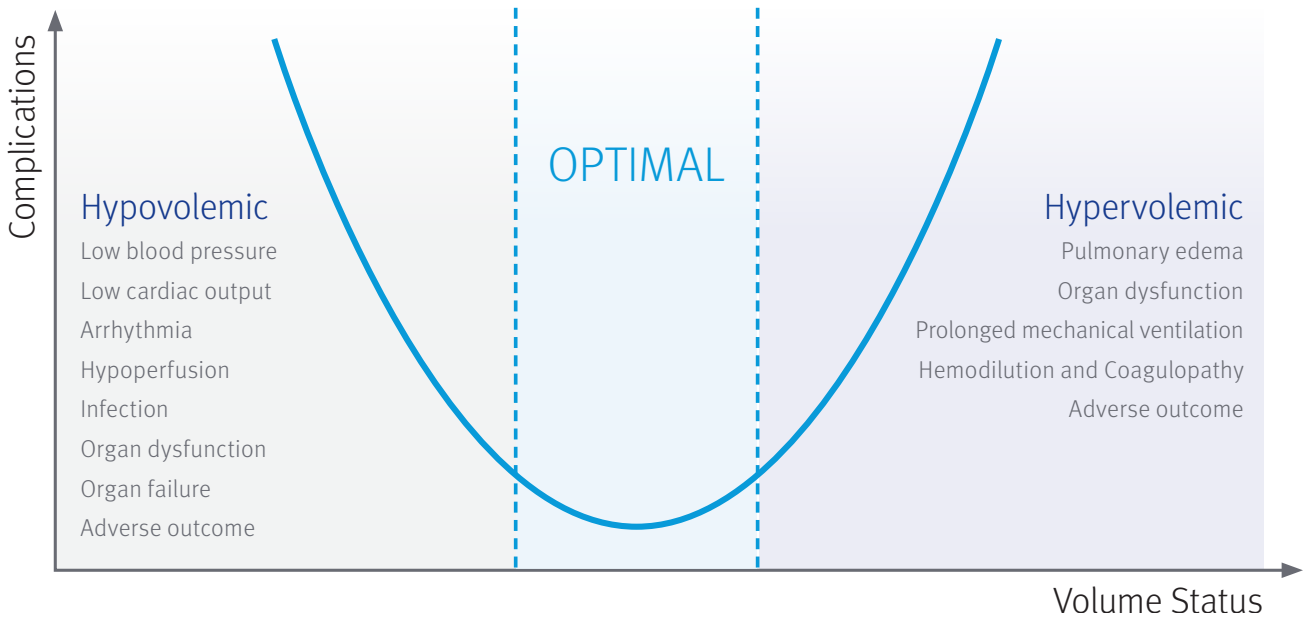

Intraoperative fluid management protocols are designed to maintain a patient’s preoperative euvoemia, or properly hydrated state by avoiding excess salt and water as well as hypovolemia conditions known to increase postoperative complication rates.^{18,19,20,21} Once optimal intravascular volume is established, intravenous fluid administration is controlled and increased only if clinically indicated.

Intravascular hypovolemia can lead to complications stemming from hypoperfusion of vital organs and the bowel. However, administering too much fluid can lead to bowel edema and increased interstitial lung water, resulting in complications and delayed return of gastrointestinal function.^{18,20,22} As a result, intravascular volume is one of the key determinants of cardiac output (CO) and oxygen delivery to the tissues. Determining the correct amount of fluid required is simplified with the focus on easily accessible, flow-based hemodynamic parameters.

To achieve optimal fluid balance, clinicians must monitor a variety of hemodynamic targets, including cardiac index (CI), stroke volume (SV), stroke volume variation (SVV), and pulse pressure variation (PPV).

MAINTAINING FLUID BALANCE WITH ACCURATE HEMODYNAMIC MONITORING

Clinicians have studied both restrictive and liberal fluid regimens as a means to guide PGDT, finding that it is important to evaluate patients on the basis of being in a state of “fluid balance.”²⁰ To achieve this optimal fluid balance, clinicians must monitor a variety of hemodynamic targets, including cardiac index (CI), stroke volume (SV), stroke volume variation (SVV), and pulse pressure variation (PPV).

Routine hemodynamic measurements, such as heart rate and mean arterial pressure (MAP), remain relatively unchanged despite reduced blood flow and are considered insensitive indicators of hypovolemia or changes in CI.^{23,24} As a result, conventional fluid management is based on clinical assessment, vital signs, central venous pressure (CVP) monitoring, or a combination of these. However, recent studies have shown that CVP is not able to predict fluid responsiveness nor can changes in blood pressure be used to approximate changes in SV or CO.^{25,26}

Many of the clinical parameters that help clinicians effectively manage patient fluids are considered dynamic due to their continuous updating and sensitivity to fluid volume changes in patients. These dynamic flow parameters, such as SV, SVV, and PPV, are needed to accurately predict fluid responsiveness.²⁷ More than two dozen randomized controlled studies and meta-analyses support the advantages of hemodynamic optimization over standard fluid management to achieve superior clinical outcomes.^{12,28,29,30}


DEMONSTRATING THE CLINICAL EFFECTIVENESS OF PGDT

PGDT has proven to help improve outcomes for a variety of patient populations, including orthopedic cases, which often carry a significant burden of complications and mortality when performed in elderly patients. One study examining the effectiveness of PGDT protocols when used with total hip replacement (THR) surgery patients found that PGDT patients received a greater volume of fluids during the intraoperative period, exhibited higher urine output, achieved a more positive fluid balance and ultimately experienced fewer postoperative complications.¹⁴ In this study, SV and CO were tracked continuously and oxygen delivery index (DO₂I) was calculated by inputting the hemoglobin concentration and SaO₂ into standard equations.

PGDT based on pulse pressure variation (PPV) has also been shown to improve patient outcomes while utilizing less invasive, pulse contour-based, arterial pressure monitoring devices to determine CO.^{31,32,33} These devices provide a wide variety of dynamic parameters, are generally easier to use, and have exhibited the ability to adequately assess dynamic changes in CO.³⁴ PPV, which is inversely proportional to SV when attempting to maximize SV intraoperatively, has also been shown to accurately reflect volume responsiveness in a number of different high-risk surgical groups, helping avoid unnecessary and potentially harmful volume loading by reflecting the cyclic changes in preload induced by mechanical ventilation.^{35,36,37,38}

A randomized, multi-center study of patients undergoing major abdominal surgery demonstrated that PGDT using PPV, radial artery pulse contour CI, and MAP leads to a reduction in postoperative complications (See Figure 2).³⁹ After an initial assessment, patients in the study were reassessed every 15 minutes intraoperatively to maintain hemodynamically stable values. Hemodynamic data were documented every 30 minutes, and ventilatory parameters every 60 minutes. Meanwhile, the treatment of patients in the control group was performed at the discretion of the care-giving anesthesiologist.

FIGURE 2. MULTI-CENTER COMPARISON OF COMPLICATIONS FOLLOWING MAJOR SURGERY BETWEEN PGDT AND NON-PGDT GROUPS³⁹

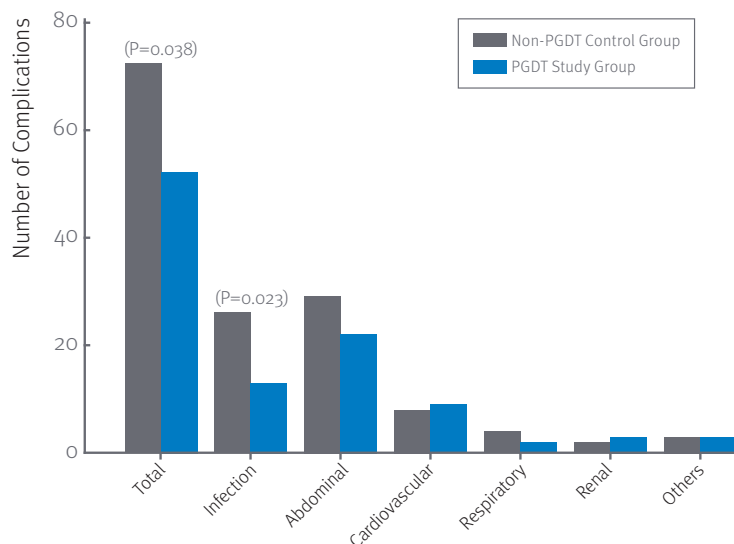

Control group (basic anesthesia monitoring without PPV and CI information; treatment at the discretion of the caregiver); Study group (basic anesthesia monitor including continuous PPV and CI information; initial hemodynamic assessment/intervention based on PPV and CI treatment algorithms)


Studies show that applying ERAS practice guidelines reduces postoperative complications by up to 50% and hospital LOS by 2.5 days.^{28,44}

Study data showed that the number of patients experiencing a complication, such as renal failure, respiratory failure, or wound infections was significantly lower in the PGDT study group.

One study examining the effects of PGDT based on an SVV-guided protocol (See Figure 3) showed that by maintaining an SVV of <12% helped reduce complications and LOS for moderate-risk patients following major surgery.⁴⁰ In the study, the PGDT group experienced earlier return of GI function (3 vs. 4 days) and oral intake (4 vs. 5 days), decreased hospital stay (5 vs. 7.5 days), and significantly higher recovery scores compared to the non-PGDT control group. Another study of patients undergoing major surgeries showed that monitoring SV increases greater than 10% with 250 mL boluses reduced total LOS by 3.6 days (See Figure 4).⁴¹


If large amounts of blood loss with severe hemodynamic instability are anticipated during the course of surgery, such as vascular surgery or liver transplant, arterial pressure-based hemodynamic monitoring may be warranted.³⁹ Research shows that, for some high-risk patients, the use of central venous oxygen saturation (ScvO₂) in the early postoperative period is an effective parameter to help guide fluid therapy.^{42,43}

FIGURE 3. GDT FLUID MANAGEMENT PROTOCOL⁴⁰

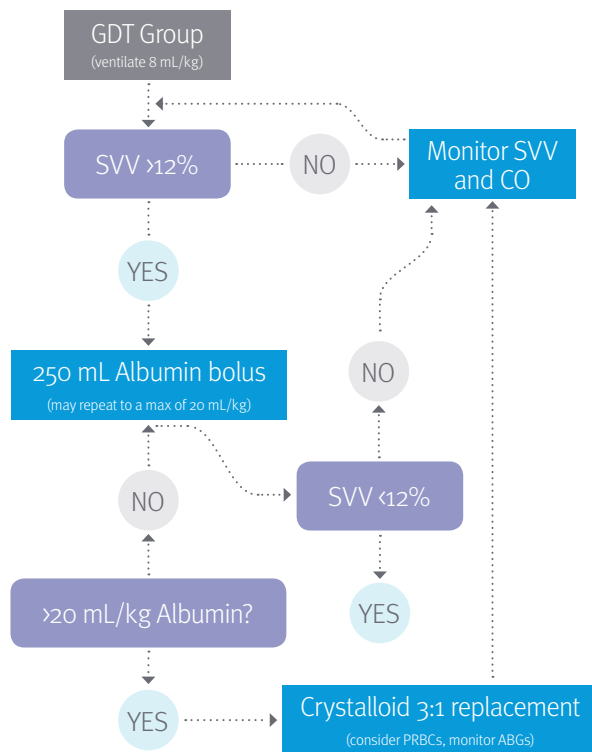

The above SVV-based GDT protocol was used to manage patients by maintaining an SVV <12%. Data from the study indicates this intraoperative protocol, guided by arterial pressure-based cardiac output monitoring, may impact postoperative outcomes and hospital LOS.

FIGURE 4. STROKE VOLUME OPTIMIZATION PROTOCOL⁴¹

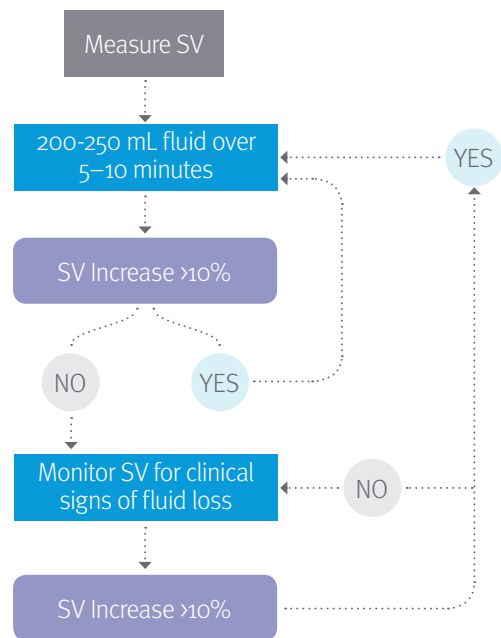

Abbreviation: SV = Stroke Volume

The above SV optimization protocol has been associated with reduced postoperative complications and hospital LOS and is now officially recommended by the National Institute for Clinical Excellence in the UK and by the French Society of Anesthesiology & Intensive Care (SFAR).

One review of twenty-nine randomized controlled trials (4,805 patients) found that proactive hemodynamic intervention in the perioperative period for cardiovascular management of moderate- to high-risk patients was associated with a 52% reduction in mortality in studies using a pulmonary artery catheter as well as those that employed cardiac index or oxygen delivery as the end-point targets.³⁰ The meta-analysis also found that preemptive hemodynamic intervention was associated with a 57% reduction in overall rates of surgical complications.

ECONOMIC BENEFITS

Meta-analyses of published studies focused on major abdominal surgery show that applying ERAS practice guidelines reduces postoperative complications by up to 50% and hospital LOS by 2.5 days.^{28,44} One study evaluating the economic impact of the clinical benefits of ERAS found that both direct medical and indirect non-medical costs were significantly lower in the ERAS group.⁴⁵ A similar published cost analysis of ERAS for colorectal surgery showed that

the full financial burden of setting up and maintaining an ERAS program was significantly offset by the costs saved by reducing postoperative resource utilization.⁴⁶ For the 50 patients managed using ERAS protocols, clinicians documented a significant reduction in total hospital stay, intravenous fluid use, complications, and duration of epidural use, representing an overall cost-savings of roughly \$345,000, or \$6900 per patient—more than offsetting the \$102,000 cost associated with implementing the ERAS program.

Although ERAS protocols were initially created for colorectal surgery, the concept has been studied and adapted for use across a wide range of medical specialties including gynecology, thoracic, vascular, pediatric, and orthopedic surgery.

According to a study conducted at the University of Virginia Medical Center, ERAS patients undergoing colorectal surgery reduced their length of hospital stay to 4.5 days compared to 6.9 days for patients who didn't experience ERAS protocols.⁴⁷ The hospital also cut its surgical complication rate from 30% to 15% and saved an average of \$7,000 per patient, totaling \$700,000 over six months. ERAS programs also allowed the hospital to accept more transfer patients with added bed capacity resulting from shorter LOS. In addition, patient engagement and satisfaction increased significantly, and the 30-day readmission rate fell from 19% to 7%, suggesting the benefits of ERAS protocols extend beyond the hospital walls.

CONCLUSION

The implementation of PGDT protocols guided by continuous hemodynamic monitoring is associated with improved postoperative outcomes and decreased LOS. ERAS guidelines were initially created for colorectal surgery but have since been studied and adapted for use across a wide range of medical specialties, including gynecology, thoracic, vascular, pediatric, and orthopedic surgery. The current evidence base shows that perioperative management, specifically the use of PGDT guided by real-time, continuous hemodynamic monitoring, helps clinicians maintain a patient's optimal fluid balance. Clinical awareness of the impact of ERAS is continuing to grow with help from organizations and societies such as the Enhanced Recovery Partnership (ERP), the Association of Surgeons of Great Britain and Ireland, the French Society of Anesthesiology (SFAR), and the Enhanced Recovery After Surgery Society.^{48,49,50,51}

References

1. Khuri SF, Henderson WG, DePalma RG, Mosca C, Healey NA, Kumbhani DJ; Participants in the VA National Surgical Quality Improvement Program. Determinants of long-term survival after major surgery and the adverse effect of postoperative complications. *Ann Surg* 2005;242:326–41
2. Schilling, et al. Prioritizing Quality Improvement in General Surgery. *J Am Coll Surg*. 2008. (n=129,233)
3. Varadhan KK, Lobo DN, Ljungqvist O. Enhanced Recovery After Surgery: The Future of Improving Surgical Care. *Crit Care Clin*. 2010 Jul;26(3):527–47
4. Eskicioglu C, Forbes SS, Aarts MA, et al. Enhanced recovery after surgery (ERAS) programs for patients having colorectal surgery: a meta-analysis of randomized trials. *J Gastrointest Surg*. 2009;13:2321–9.
5. Basse L, Raskov HH, Hjort Jakobsen D, Sonne E, Billesbolle P, Hengel HW, et al. Accelerated postoperative recovery programme after colonic resection improves physical performance, pulmonary function and body composition. *Br J Surg* 2002;89(4):446–53.
6. Muller S, Zalunardo MP, Hubner M, Clavien PA, Demartines N. A fast-track program reduces complications and length of hospital stay after open colonic surgery. *Gastroenterology* 2009;136(3):842–7.
7. Miller T, Thacker J, White W, et al. Reduced length of hospital stay in colorectal surgery after implementation of an enhanced recovery protocol. *Anesth Analg*. 2014 May;118(5):1052–61.
8. Podore PC, Throop EB. Infra-renal aortic surgery with a 3-day hospital stay: A report on success with a clinical pathway. *J Vasc Surg*. 1999;29:787–92.
9. Tovar EA, Roethe RA, Weissig MD, et al. One-day admission for lung lobectomy: an incidental result of a clinical pathway. *Ann Thorac Surg*. 1998;65:803–6.
10. Koupparis A, Dunn J, Gillatt D, et al. Improvement of an enhanced recovery protocol for radical cystectomy. *British Journal of Medical and Surgical Urology*. 2010;3:237–40.
11. Cecconi M, Fasano N, Langiano N, Divella M, Costa MG, Rhodes A, Della Rocca G. Goal-directed hemodynamic therapy during elective total hip arthroplasty under regional anesthesia. *Crit Care* 2011, 15:R132.
12. Giglio MT, Marucci M, Testini M, Brienza N. Goal-directed haemodynamic therapy and gastrointestinal complications in major surgery: a meta-analysis of randomized controlled trials. *British Journal of Anaesthesia* 2009 Nov;103(5):637–46.
13. Bellamy, et al. Wet, Dry, or Something Else? Editorial. *BJA*, 2006.
14. Sinclair S, James S, Singer M: Intraoperative intravascular volume optimization and length of hospital stay after repair of proximal femoral fracture: randomized controlled trial. *BMJ* 1997, 315:909–912.
15. Bartha E, Arfwedson C, Imnell A, Fernlund ME, Andersson LE, Kalman S. Randomized controlled trial of goal-directed hemodynamic treatment in patients with proximal femoral fracture. *Br J Anaesth* 2013, 110:545–553.
16. Noblett SE, Snowden CP, Shenton BK, Horgan AF. Randomized clinical trial assessing the effect of Doppler-optimized fluid management on outcome after elective colorectal resection. *Br J Surg* 2006, 93:1069–1076.
17. Cannesson M, Ramsingh D, Rinehart J, et al. Perioperative goal-directed therapy and postoperative outcomes in patients undergoing high-risk abdominal surgery: a historical-prospective, comparative effectiveness study. *Critical Care* 2015, 19:261.
18. Brandstrup B, Tonnesen H, Beier-Holgersen R, Hjortso E, Ording H, Lindorff-Larsen K et al. Effects of intravenous fluid restriction on postoperative complications: comparison of two perioperative fluid regimens: a randomized assessor-blinded multicenter trial. *Ann Surg* 2003; 238: 641–648.
19. Lobo DN, Bostock KA, Neal KR, Perkins AC, Rowlands BJ, Allison SP. Effect of salt and water balance on recovery of gastrointestinal function after elective colonic resection: a randomised controlled trial. *Lancet* 2002; 359:1812–1818.
20. Varadhan KK, Lobo DN. A meta-analysis of randomized controlled trials of intravenous fluid therapy in major elective open abdominal surgery: getting the balance right. *Proc Nutr Soc* 2010; 69: 488–498.
21. Strunden MS, Heckel K, Goetz AE, Reuter DA. Perioperative fluid and volume management: physiological basis, tools and strategies. *Ann Intensive Care* 2011, 1(1):2.
22. Chowdhury AH, Lobo DN. Fluids and gastrointestinal function. *Curr Opin Clin Nutr Metab Care* 2011;14:469–76.
23. Mythen MG, Webb AR: Intraoperative gut mucosal hypoperfusion is associated with increased post-operative complications and cost. *Intensive Care Med* 1994, 20:99–104.
24. Pierrakos C, Velissaris D, Scolletta S, Heenen S, De Backer D, Vincent J-L. Can changes in arterial pressure be used to detect changes in cardiac index during fluid challenge in patients with septic shock? *Intensive Care Med* 2012, 38:422–428.
25. Marik PE, Baram M, Vahid B. Does central venous pressure predict fluid responsiveness? A systematic review of the literature and the tale of seven mares. *Chest*. 2008 Jul;134(1):172–8.
26. Le Manach Y, Hofer CK, Lehot JJ, Vallet B, Goarin JP, Tavernier B, Cannesson M. Can changes in arterial pressure be used to detect changes in cardiac output during volume expansion in the perioperative period? *Anesthesiology*. 2012 Dec;117(6):1165–74.
27. Michard F, Biais M. Rational fluid management: dissecting facts from fiction. *Br J Anaesth*. 2012 Mar;108(3):369–71.
28. Corcoran T, Rhodes JE, Clarke S, Myles PS, and Ho KM. Perioperative fluid management strategies in major surgery: a stratified meta-analysis. *Anesth Analg*. 2012 Mar;114(3):640–51.
29. Dalfino L, Giglio MT, Puntillo F, Marucci M, Brienza N: Hemodynamic goal-directed therapy and postoperative infections: earlier is better. A systematic review and meta-analysis. *Crit Care* 2011, 15:R154.
30. Hamilton MA, Cecconi M, and Rhodes A. A systematic review and meta-analysis on the use of preemptive hemodynamic intervention to improve postoperative outcomes in moderate and high-risk surgical patients. *Anesth Analg*. 2011 Jun;112(6):1392–402.
31. Benes J, Chytra I, Altmann P, Hluchy M, Kasal E, Svitak R, Pradl R, Stepan M. Intraoperative fluid optimization using stroke volume variation in high risk surgical patients. *Crit Care* 2010, 14:R118.
32. Lopes MR, Oliveira MA, Pereira VO, Lemos IP, Auler JO Jr, Michard F. Goal-directed fluid management based on pulse pressure variation monitoring during high-risk surgery: a pilot randomized controlled trial. *Crit Care* 2007, 11:R100.
33. Breukers RM, Seppehrkhouy S, Spiegelenberg SR, Groeneveld AB. Cardiac output measured by a new arterial pressure waveform analysis method without calibration compared with thermodilution after cardiac surgery. *J Cardiothorac Vasc Anesth* 2007, 21:632–635.
34. Petzoldt M, Riedel C, Braeunig J, Haas S, Goepfert MS, Treede H, Baldus S, Goetz AE, Reuter DA. Stroke volume determination using transcatheter pulmonary thermodilution and arterial pulse contour analysis in severe aortic valve disease. *Intensive Care Med* 2013, 39:601–611.
35. Biasi M, Bernard O, Ha JC, Degryse C, Sztark F. Abilities of pulse pressure variations and stroke volume variations to predict fluid responsiveness in prone position during scoliosis surgery. *Br J Anaesth* 2010, 104:407–413.
36. Derichard A, Robin E, Tavernier B, Costecalde M, Fleyfel M, Onimus J, Lebuffe G, Chambon JP, Vallet B. Automated pulse pressure and stroke volume variations from radial artery: evaluation during major abdominal surgery. *Br J Anaesth* 2009, 103:678–684.
37. Hadian M, Severyn DA, Pinsky MR. The effects of vasoactive drugs on pulse pressure and stroke volume variation in postoperative ventilated patients. *J Crit Care* 2011, 26:328.e1–8.
38. Reuter DA, Felbinger TW, Kilger E, Schmidt C, Lamm P, Goetz AE. Optimizing fluid therapy in mechanically ventilated patients after cardiac surgery by on-line monitoring of left ventricular stroke volume variations. Comparison with aortic systolic pressure variations. *Br J Anaesth* 2002, 88:124–126.
39. Salzwedel C, Puig J, Carstens A, Bein B, Molnar Z, Kiss K, Hussain A, Belda J, Kirov MY, Sakka SG, Reuter DA. Perioperative goal-directed hemodynamic therapy based on radial arterial pulse pressure variation and continuous cardiac index trending reduces postoperative complications after major abdominal surgery: a multi-center, prospective, randomized study. *Crit Care*. 2013 Sep 8;17(5):R191.
40. Ramsingh, D.S., Sanghvi, C., Gamboa, J, et al, Outcome impact of goal directed fluid therapy during high risk abdominal surgery in low to moderate risk patients: a randomized controlled trial. *J Clin Monit Comput*.2013;27:249–257.
41. Kuper M, Gold SJ, Callow C, et al. Intraoperative fluid management guided by oesophageal Doppler monitoring. *BMJ*. 2011;342:d3016.
42. Futier E, Constantin JM, Petit A, Chanques G, Kwiatkowski F, Flamein R, et al. Conservative vs restrictive individualized goal-directed fluid replacement strategy in major abdominal surgery: a prospective randomized trial. *Arch Surg* 2010;145(12):1193–200.
43. Pearse R, Dawson D, Fawcett J, Rhodes A, Grounds RM, Bennett ED. Changes in central venous saturation after major surgery, and association with outcome. *Crit Care* 2005;9(6):R694–9
44. Varadhan KK, Neal KR, Dejong CH, Fearon KC, Ljungqvist O, Lobo DN. The enhanced recovery after surgery (ERAS) pathway for patients undergoing major elective open colorectal surgery: a meta-analysis of randomized controlled trials. *Clinical Nutrition (Edinburgh, Scotland)* 2010 Aug;29(4):434–40.
45. King PM, Blazeby JM, Ewings P, et al. Randomized clinical trial comparing laparoscopic and open surgery for colorectal cancer within an enhanced recovery programme. *Br J Surg*. 2006;93:300–8.
46. Sammour T, Zargar-Shoshtari K, Bhat A, et al. A programme of Enhanced Recovery After Surgery (ERAS) is a cost-effective intervention in elective colonic surgery. *NZ Med J*. 2010 Jul 30;123(1319):61–70.
47. Thiele RH, Rea KM, Turrentine FE, Friel CM, Hassinger TE, Goudreau BJ, Umapathi BA, Kron IL, Sawyer RG, Hedrick TL, McMurry TL. Standardization of care: impact of an enhanced recovery protocol on length of stay, complications, and direct costs after colorectal surgery. *J Am Coll Surg*. 2015 Apr;220(4):430–43.
48. Mythen et al. Perioperative fluid management: Consensus statement from the enhanced recovery partnership. *Perioperative Medicine* 2012.
49. Powell-Tuck, et al. British consensus guidelines on intravenous fluid therapy for adult surgical patients.
50. Vallet, et al. Guidelines for perioperative haemodynamic optimization. *Ann Fr Anesth Reanim* 2013.
51. Gustafsson, et al. Guidelines for perioperative care in elective colonic surgery: enhanced recovery after surgery (ERAS) society. *World J Surg* 2013.

